



A 100% Australian owned company supplying scientific laboratories throughout Australia and overseas since 1987.



Just some of the latest and most popular laboratory products available from Rowe Scientific Pty Ltd.

Timer Stopwatch /clock IT0240

- Displays hours, minutes and seconds
- Stopwatch function
- Count-down timer counts up once alarm goes off to show overtime.
- Clock can be set as either 12 hour am/pm or 24 hour format
- 24 hour timer capacity with 1 second resolution
- Timer memory allows a preset time to be stored
- Magnetic back clip
 - a) adheres to metal surfaces
 - b) clips on to pockets and the like, or
 - c) operates as "easel" type stand
- Stopwatch, count-down and clock can operate simultaneously
- Battery included



Scan code to order.



List \$27.00 now
\$18.90

Call for bulk buy discount

Thermometer max/min inside/outside IT0248

- Max/Min memory records high and low temperature
- Measures inside and outside temperature through the use of two sensors
- Wide temperature range of -50°C to +70°C
- Temperature readout in 0.1°C
- Three metre long weather-proof sensor wire for outside temperature readings
- Includes wall mount and desktop stand



Scan code to order.



List \$40.50 now
\$24.30

Call for bulk buy discount



- ✓ One Stop Shop
- ✓ 100% Australian
- ✓ Online ordering
- ✓ Phone ordering
- ✓ Account Managers
- ✓ Leading Brands
- ✓ Chemicals
- ✓ Consumables
- ✓ Equipment

Below are just some of the leading brands we offer.



EX-DEMO MODEL - BE QUICK!

Mettler Toledo Precision balance, New Classic ML3002T with touchscreen 3200 x 0.01g - **1 ONLY**
Rowe Code IB0174



INTERNATIONAL YEAR OF THE PERIODIC TABLE www.iypt2019.org

SEASON GREETINGS





PARAFILM[®] M Cutter

Rowe code: CB2023




Easy to set-up
Easy to use
Efficient and time saving

1. Before opening the cutter move the blade carriage to the stop on the left. Lift both lateral tongues to open.

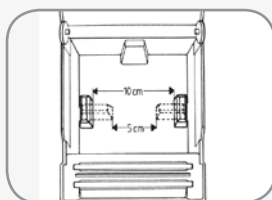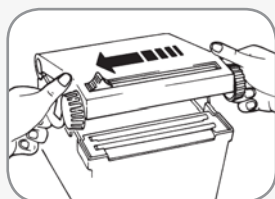

2. Adjust the fitting brackets to the width of the roll.

3. Open the film carton and pull out the film. Insert the box into the cutter.

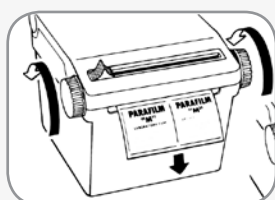

4. Turn the handwheels backwards until desired length of film appears. Cut the film by moving the blade carriage across.

List \$426.00 now

\$355.00

and don't forget the film ...

PP1200
Parafilm 5cm wide x 75m long

PP1250
Parafilm 10cm wide x 38m long



Filter Paper



FILTER PAPER 11.0cm similar to Gr41 (MN640w)
FP0037 Unit: 100 List \$14.36 **SALE \$12.21**

FILTER PAPER 11.0cm similar to Gr1 (MN615)
FF1200 Unit: 100 List \$16.24 **SALE \$13.80**

FILTER PAPER 11.0cm similar to Gr40 (MN 640 md)
FF2524 Unit: 100 List \$54.47 **SALE \$46.30**

FILTER PAPER 12.5cm similar to Gr4 (MN617)
FF2220 Unit: 100 List \$25.33 **SALE \$21.53**

FILTER PAPER 12.5cm similar to Gr541 (MN1640W)
FF3600 Unit: 100 List \$66.11 **SALE \$56.19**

FILTER PAPER 15.0 cm similar to Gr91 (MN126/70)
FP0052 Unit: 100 List \$14.55 **SALE \$12.37**

FILTER PAPER 15.0cm similar to Gr91 (MN751)
CREPED **FF3050** Unit: 100 List \$21.01 **SALE \$17.86**

GLASS FIBRE FILTER PAPER

FILTER PAPER GLASS FIBRE 4.7cm similar to Gr GF/A (MNGF-1) **FF4627** Unit: 100 List \$29.48 **SALE \$25.06**

FILTER PAPER GLASS FIBRE 4.7cm similar to Gr GF/C (MN GF-3) **FF4128** Unit: 100 List \$36.29 **SALE \$30.85**

FILTER PAPER GLASS FIBRE 5.5cm SIMILAR Gr GF/C (MN GF-3) **FF4144** Unit: 100 List \$45.36 **SALE \$38.56**

FILTER PAPER GLASS FIBRE 7.0cm 0.6micron (MN GF/3) **FF4150** Unit: 100 \$67.66 **SALE \$57.51**

FILTER PAPER GLASS FIBRE 9.0cm GF-1 **FF4601** Unit: 100 List \$68.80 **SALE \$58.48**

Macherey-Nagel are a German manufacturer of quality filter papers and more since 1911. Above is only a sample of what we hold. Call us to discuss your needs and special pricing on bulk orders.





A 100% Australian owned
company since 1987.

OCTOBER - DECEMBER 2019

ROWE SCIENTIFIC PTY LTD

LABLIKES

Are you using the right gloves?

**NeoTouch® 25-101 / 201 Neoprene
POWDER-FREE, TEXTURED
Unique Glove combining Allergy Prevention
and Comfort.**

Excellent resistance to acids, bases and alcohols. Proprietary formulation provides superior comfort. Excellent donning due to special coating. Secure grip in both dry and wet conditions. Thickness: 0.13mm. Length: 240 or 285mm. 10 X 100 in a box.

**Microflex® 93-852 Nitrile
POWDER-FREE, DISTINCT BLACK COLOUR.
Strong and versatile**

Colour masks stains and provides a contrasting backdrop for lighter material. Exceptional flexibility with high tensile strength providing lasting protection. Fully textured for a consistent firm grip. Non-foaming formula allows for exceptional wet grip. Made without natural rubber latex to help Type 1 skin allergies, skin irritation and prevent dryness. Thickness: 0.12mm Length:240mm. 100 per box.

**Microflex® 93-260 Nitrile & Neoprene
TOUGH CHEMICAL PROTECTION,
UNPARALLELED COMFORT. Thinnest
disposable chemical resistant glove**

Triple layer design for protection against harsh chemicals (eg heptane, formaldehyde, sulfuric acid). Thin construction that enhances tactility and dexterity. Soft polymer formulation and design that improves fit, feel and flexibility for more comfort and extended wear time. Low acceptable pinhole rate (0.65 AQL) and extended cuff for reliable protection against hazardous substances. Thickness: 0.198mm, Length: 300mm. 50 in a box.

**TouchNTuff® 92-600 / 605 Nitrile
POWDER-FREE, SMOOTH / TEXTURED
Best in Class Splash Resistant Nitrile Glove.**

Easy donning and strong grip in wet or dry conditions. High comfort and dexterity -Thin nitrile technology. Outstanding protection against chemical splashes. Robust glove – excellent puncture resistance. Thickness: 0.12mm Length: 240mm. 100 in a box.



**NeoTouch® 25-201
HG0240**

**Microflex® 93-852
HG0081**

**MicroFlex® 93-843
HG0190**

**Microflex® 93-260
HG0162**

**TouchNTuff® 92-600
HG1897**

**Microflex® 93-833
HG0153**

**MicroFlex® 93-843/853 Nitrile.
POWDER-FREE, TEXTURED. Sturdy glove
for demanding applications**

Provides durable, heavy-duty coverage. Textured fingertips. Lowest pinhole rates (AQL = 0.65). Thickness: 0.11mm. Length: 245mm. 100 in a box.

**Microflex® 93-833 Nitrile
POWDER-FREE, ERGONOMIC DESIGN. Thinnest
and strongest ergonomically certified glove**

Textured fingertips. Lowest pinhole rates (AQL 0.65). Reduces hand fatigue and muscle effort (US ergonomics certification). Thickness: 0.07mm. Length: 240mm. 250 in a box.

AND THERE'S MORE.
Speak with your local Rowe Scientific consultant
to review your glove requirements.



Most gloves available in
sizes ranging from:
XS, S, M, L, XL, XXL

Need to shake things up?



IV1085

MIXER VORTEX ZX3 CLASSIC with touch sensor or continuous mode. 0-3000 rpm



List \$354.49 now

\$272

IM1081

MIXER, VORTEX CLASSIC with IR sensor, TX4. 0-3000 rpm, with display.



List \$599.72 now

\$462

IM1063

MIXER VORTEX, WIZARD with IR sensor tube detection or continuous mode. Low profile. 0-3000



List \$492.88 now

\$383

Micro Tubes

VOL	DESCRIPTION	UNIT	PRICE	ROWE CODE
15mL	STERILE conical bottom 120mmH x 16mmD Polypropylene	500	95.00	PT1771
50mL	STERILE conical bottom 115mmH x 28mmD Polypropylene	500	115.00	PT1170
15mL	NON-STERILE conical bottom 17mmD x 120mmH Graduated Polypropylene with grey cap	500	78.75	PT0444
50mL	NON-STERILE conical bottom 115mmH x 30mmD Graduated Polypropylene with blue cap. Tubes and caps separate	500	90.00	PT0099



IM1064 Mixer vortex accessory, foam stand for 19 x 1.5mL Eppendorf tubes suits Velp Classic and Wizard.



IM1065 Mixer vortex accessory, foam stand soft top customisable suits Velp Classic and Wizard




IM1066 Mixer vortex accessory, foam stand for 5x16mm test tubes suits Velp Classic and Wizard



IM1072 Mixer vortex accessory, foam stand for 4x29mm test tubes



IM1067 Mixer vortex accessory, foam stand for microtitre plate suits Velp Classic and Wizard

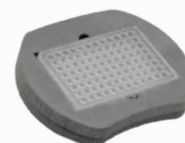

IM1068 Mixer vortex accessory, rubber supporting plate 50mm diameter suits Velp Classic and Wizard



IV1074 Black rubber cup to suit Wizard-Classical-WX vortex mixer (new improved rubber compound)





SPL LIFE SCIENCES PRODUCTS



SPL Culture Flasks

SPL culture flasks are ergonomically designed for easy handling, with minimal chance of contamination. Short, wide and angled neck design enables easy access of serological pipettes and cell scrapers.

- Ergonomic design to facilitate easy handling and minimise contamination
- Short & wide neck with angled design allows easy access
- Excellent stackability
- Plugs & filter caps are available for all flasks
- Available in a range of sizes: 25cm², 75cm², 175cm²
- Venting position for 175cm² plug cap flasks
- Non-treated flasks (for suspended cultures) are provided with white caps
- Non-pyrogenic
- Non-cytotoxic
- DNase/RNase-free
- Human DNA-free



SURFACE AREA	CELL CULTURE FLASKS	PACK SIZE	PACK PRICE	ROWE CODE
25cm ²	T25 SPL POLYSTYRENE with HPDE filter cap (vented) TC treated sterile. 40 bags x 5 flasks	200	129.59	PC0329
75cm ²	T75 SPL POLYSTYRENE with HPDE filter cap (vented) TC treated sterile. 20 bags x 5 flasks	100	125.69	PC0330
175cm ²	T175 SPL POLYSTYRENE with HPDE filter cap (vented) TC treated sterile. 8 bags x 5 flasks	40	100.55	PC0331

NOTE: Bulk buy discounts are available for volume users.



SPL Serological pipettes (Suitable for cell culture work)

SPL Serological pipettes are classified by sample volume. Ascending & descending graduations facilitate the reading of both dispensing and remaining volume.

- Six different capacities/sizes
- Colour-coded for easy identification
- Various packaging sizes
- Non-pyrogenic
- Non-cytotoxic
- DNase/RNase-free
- Human DNA-free

SEROLOGICAL PIPETTES	PACK SIZE	PACK PRICE	ROWE CODE
1mL SPL Yellow Individual Sterile	800	158.01	PS0332
2mL SPL Green Individual Sterile	600	132.87	PS0333
5mL SPL Blue Individual Sterile	400	112.52	PS0334
10mL SPL Orange Individual Sterile	400	122.10	PS0335
25mL SPL Red Individual Sterile	200	92.78	PS0336
50mL SPL Purple Individual Sterile	160	287.28	PS0337

This is just a sample of the SPL range of life sciences products available from Rowe Scientific PL. If you are seeking something specific please contact your local Rowe consultant.



Rowe pH Buffers

Rowe Scientific introduced the unique twin neck bottle design for pH buffers to the Australian scientific community.

Features include:

- Separate sampling chamber reduces contamination.
- Small volume can be discarded once complete
- Overuse is reduced
- As used in European laboratories

All our buffers are traceable to NBS (NIST) international standards.



Royal Flying Doctor Service

In support of the Royal Flying Doctor Service of Australia, every purchase of Rowe pH buffers (CB3515, CB3115, CB2660) in the unique twin neck bottles, results in a \$2.00 per item donation. Scan the QR code to see our RFDS meter.



CB3115

Buffer pH 7.00 (Green) solution NIST traceable, in contamination proof bottle.

500mL

1x \$20.41 ea

6x \$11.41 ea

CB2660

Buffer pH 4.00 (Red) solution NIST traceable, in contamination proof bottle.

500mL

1x \$20.75 ea

6x \$11.23 ea

CB3515

Buffer pH 10.00 (Blue) solution NIST traceable, in contamination proof bottle.

500mL

1x \$19.68 ea

6x \$10.82 ea

Buffer also available in larger 1L, 5L and 20L containers.



Centrifuge Deals

Buy a centrifuge below and receive a **FREE COMPATIBLE ROTOR VALUED AT OVER \$900**



IC0008

CENTURION K1015 Microtube high speed centrifuge

Normally \$3348

SPECIAL PRICE \$2350

SAVE \$998

- Accepts 0.2 - 2.2mL micro tubes
- 500 - 15,000 rpm (1 rpm steps)
- RCF max of 22,000 x g
- Powerful 172 watt motor
- 0 - 99 minute timer (1 second steps)
- 235 x 235 x 350mm dimensions
- 10 inbuilt acceleration and deceleration programs
- 108 program memory



Includes free IR4110 Rotor valued at \$900



IC0103

CENTURION C1015 Pro-analytical high speed centrifuge

Normally \$2471

SPECIAL PRICE \$1990

SAVE \$481

- Accepts 0.2 - 15mL tubes
- 500 - 15,000 rpm (10 rpm steps)
- RCF max of 22,000 x g
- Powerful 140 watt motor
- 0 - 99 minute timer (30 second steps)
- 235 x 235 x 350mm dimensions
- 10 program memory



Includes free IR4136 Rotor valued at \$480

Stocks are limited so be quick.

Offer available while stocks last.



Contact us today.



Laboratory Sterilisers



LABOKLAV 25

Compact benchtop design with unique viewing window. Laboklav offer a range of solutions. Speak to one of our consultants about your steriliser needs.

IL0086 Autoclave bench top Laboklav 25-MV air cooling with integrated vacuum pump rapid recooling steam generator automatic closing system 25L nominal capacity - can fit up to 5x 1L lab bottles 240V/10 AMP

IA0160 Autoclave bench-top basic unit Laboklav 25B with integrated steam generator automatic closing system 25L nominal capacity - can fit upto 5x 1L lab bottles 10 AMP

IA0161 Autoclave bench-top Laboklav 25M air cooling with integrated steam generator and rapid recooling automatic closing system 25L nominal capacity-can fit up to 5x 1L lab bottles 10 AMP

IA0162 Autoclave bench-top Laboklav 25V air cooling with integrated steam generator vacuum system automatic closing system 25L nominal capacity can fit up to 5x 1L lab bottles 10 AMP



LABOKLAV 55 - 195

A well designed unit at a low height for easy loading. Baskets are adapted to the chamber size to permit optimal utilisation of the chamber volume. Optional loading system and lifting devices make loading and unloading of the autoclave easy, even for heavy items requiring sterilisation.

IA0185 Autoclave Laboklav 80 MS double jacket water cooling and compressed air counter pressure 80L capacity 410mm diameter chamber x 610mmL with touch display and USB port, 3-phase power.

SHP Steriltechnik autoclaves offer a wide range of chamber sizes between 25L and 195L, universal installation options such as free-standing or tabletop units, with numerous configurations providing individual solutions for your specific requirements.

Typical areas of application for our units include all microbiological laboratories in research institutions, universities, quality control in the pharmaceutical and food industries as well as waste sterilisation, for example in the field of medical diagnosis.

Laboklav series autoclaves are characterised by their reliability and quality. Manufactured exclusively in Germany to high and continuously monitored quality standards. Units are TÜV-certified and approved for series production and bear the CE mark.

They comply with the following standards:

- Pressure Equipment Directive 97/23/EU
- DIN 58950: Sterilisation - Steam sterilisers for pharmaceutical products
- DIN 58951: Sterilisation - Steam sterilisers for laboratory use



Can also be vertically floor mounted.

BRAND™

It's a
 BRAND™
new day.



Explore the
range in our
Web Store.

Which BRAND™ dispenser do you need?

As a BRAND™ distributor we have available, and in many cases stock, a full range of dispensers to meet your laboratory requirements.

The Dispensette S and S Organic range offer both analog and digital options, with or without internal recirculation.

The Dispensette S allows for dispensing aggressive reagents directly from the bottle. Use with concentrated bases and acids such as H_3PO_4 , H_2SO_4 , HCl (<20%) and HNO_3 (<30%).

The slightly more expensive S Organic has a wider range of reagent applications. Please download and review the Chemical Compatibility Chart in the brochure link below.



Note, for Hydrofluoric acid, we recommend the use of the Dispensette S Trace Analysis bottle-top dispenser. Download and review the brochure link below.

The analogue instrument has been redesigned with an interior scalloped track which enhances setting reproducibility.



The slightly more expensive digital model allows for dynamic indication of volume setting and for calibration and adjustment according to ISO 9001 and GLP to be achieved in seconds.

All models are available with a recirculation option. This allows the dispenser to be primed after standing to ensure that the concentration of diluent is consistent, with only the minimal tip volume requiring disposal.



It's a BRAND™ new day at Rowe Scientific Pty Ltd.




A 100% Australian owned
company since 1987.

OCTOBER - DECEMBER 2019

ROWE SCIENTIFIC PTY LTD

LABLIKES

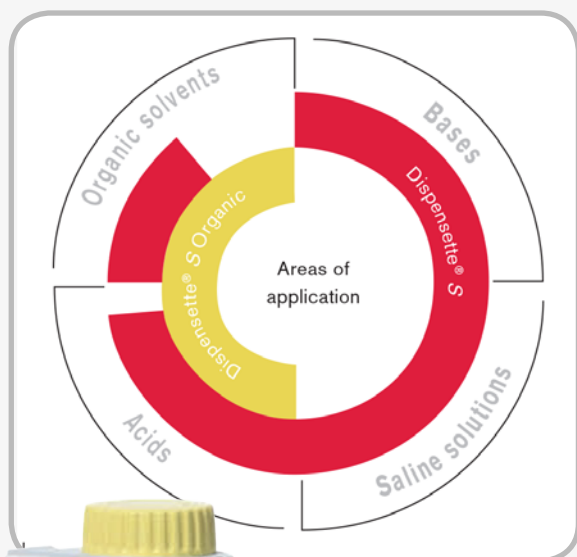

Do you require calibrated glassware?

As a distributor of BRAND™ liquid handling equipment, we are able to provide calibration certificates for your liquid handling products.



Since the beginning of 2010, DAKKS (Deutsche Akkreditierungsstelle GmbH) has been the sole provider of accreditation services for regulatory matters in Germany.

The former DKD calibration laboratory at BRAND for volumetric instruments was changed into an accredited calibration

laboratory by the DAKKS according to DIN EN ISO/IEC 17 025. DAKKS is a signatory to the International Laboratory Accreditation Cooperation (ILAC) for the mutual recognition of calibration certificates.

Please contact us if you require a calibration certificate for your new single or multi-channel pipette, burette or dispenser, glass AND plastic volumetric labware such as

- Volumetric flasks
- Pipettes
- Graduated cylinders
- Oxygen flasks (Winkler pattern)
- Density bottles

We will provide you a quote and look to deliver within 4 weeks of order.

Plan now before your next audit occurs.

We are also able to provide USP certified glassware. Selected BLAUBRAND volumetric instruments can be delivered with volume error limits in compliance with the United States Pharmacopeia (USP) with a batch certificate, or upon request, an individual certificate. The measuring instrument shows the requested error limit and the USP sign.

Rowe Scientific, helping you meet your accreditation requirements.



Submersible Stirrers



- IS0307 Stirrer magnetic Labomag 1x1 (single) position industrial standard brushed stainless steel (submersible in water up to 60°C IP68) includes separate control unit.
- IS0308 Stirrer magnetic Labomag 2x2 (4 position) with 130mm stirring point distance industrial standard brushed stainless steel (submersible in water up to 60°C IP68) includes separate control unit.
- IS0309 Stirrer magnetic Labomag 1x5 (5 position) with 65mm stirring point distance industrial standard brushed stainless steel (submersible in water up to 60°C IP68) includes separate control unit.
- IS0311 Stirrer magnetic Labomag 2x4 (8 position) with 130mm stirring point distance industrial standard brushed stainless steel (submersible in water up to 60°C IP68) includes separate control unit.



User-friendly control unit is easy to use. Integrated LEDs visualise stirring speed. The controller memorises the last setting of stirring speed and stirring direction. A USB interface is integrated as serial function.

DESCRIPTION	VOL mL	ROWE CODE	UNIT LIST PRICE EA	PACK QTY	PACK QTY SALE PRICE
Beaker low form spouted graduated borosilicate glass	25	GB1560	3.15	10	26.70
	50	GB1660	3.39	10	28.80
	100	GB1810	3.06	10	26.00
	150	GB1870	3.39	10	28.80
	250	GB2060	3.39	10	28.80
	400	GB2210	5.83	10	49.50
	600	GB2310	4.83	10	41.00
	1000	GB2610	5.74	10	48.80
	2000	GB2810	8.60	4	29.20
	3000	GB2960	57.84	4	196.60
5000	GB3110	80.11	2	136.10	
Beaker tall form spouted graduated borosilicate glass	50	GB1366	3.25	10	27.60
	100	GB3210	4.83	10	41.00
	150	GB3310	6.21	10	52.70
	250	GB3460	3.68	10	31.20
	400	GB3560	4.06	10	34.50
	600	GB3850	4.30	10	36.50
Flask erlenmeyer narrow neck rimmed borosilicate glass	1000	GB4085	5.98	10	50.80
	50	GF0227	7.83	10	66.50
	100	GF1710	3.68	10	31.20
	200	GF1826	4.83	10	41.00
	250	GF1855	4.59	10	39.00
	300	GF1986	6.07	10	51.60
	500	GF2010	6.69	10	56.80
	1000	GF2160	8.99	10	76.40
Flask erlenmeyer wide neck rimmed borosilicate glass	2000	GF2215	11.66	6	59.40
	250	GF2805	4.06	10	34.50
	300	GF2866	6.07	10	51.60
	500	GF2910	5.07	10	43.10
	1000	GF2171	9.27	10	78.80
	2000	GF2310	11.76	6	59.60

Other ranges online which may interest you ...



TPX Measuring Cylinders



Microfiltration accessories



LeadFluid Peristaltic Pumps

Or explore at www.rowe.com.au and Click the downloads tab.



How to order online at your convenience using the Rowe Scientific Web Store and its easy to use features.

As a guest, when you arrive at the Rowe Scientific Web Store, you can search for products and populate a trolley. To obtain pricing, please request a quote on your trolley of products. Clients with a login account will see pricing directly, and will be able to place orders over the Rowe Scientific Web Store.

How to Log In

If you have a login account, enter your username and password in the login panel at the top right hand corner. This will log you into your personalised Web Store account and show your special pricing.

Your Welcome Panel

When you have successfully logged in, you will see your welcome panel to the right. Here you will find some useful features that are explained on this page.

My Backorders

My backorders is a searchable list of your current back orders with Rowe Scientific. This data is updated overnight to keep you fully informed.


Quick Order

YOUR TEN MOST PURCHASED ONLINE ITEMS
Enter quantity on any of the items in the below lists then hit the 'Add all to Trolley' button at the bottom of the page.

Items	Qty	UOM
CB1525 BUFFER pH 7.0-8.0 (BULK) SOLUTION H2O2 TRACEABLE, IN CONTAMINATION PROOF BOTTLE (R608+)	579.71	500mL
FP1042 FILTER PAPER 25.0cm NET STRENGTHENED (MNT2) IN BROWN PAPER PACKING PARTICLE RETENTION IS RETENFN 7-13 MICRON (MNT33325)	548.82	100

Quick Order allows you place orders from a list of routinely ordered products. There are two parts to your quick order page. Your top ten recent purchases and a manual list made up of items you have added to your quick order list. These are grouped by product type e.g. Glassware for glass beakers etc.

To add a product to your quick order list, click on the add to quicklist button on the product listing. If the button does not appear on an item, it is already in your quick order list.

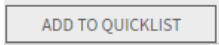

To order items directly by quick order list, simply enter the quantity required for each item being ordered then select the add to trolley button. You will see these items appear in your 'trolley' panel and the total value of your order.

You can continue searching and shopping for products or view your trolley and process your order.

Please note: only orders placed online appear in the list below

[view last 10 orders](#)

Order No. 211611 / Customer Ref. Caritest

Product ID	Description
FF2227	FILTER PAPER 15.0cm SIMILAR TO Gr4 (MN617) (434)
FF4648	FILTER PAPER GLASS FIBRE 11.0cm SIMILAR Gr G/F B

Order History

Here you can see a list of recent online orders placed via our Web Store.

Search all stock...

The Search Bar

The search bar functionality on the right allows you to string key words together as you would a Google search. e.g. 'beaker 250mL' the results will list those items found which best match the exact criteria to the top of the results page. When you use the product filter buttons a second inner search bar appears in the search results page allowing you to drill down from the results using key words.

Search by YOUR part number

We can upload your part numbers (where supplied by you) for contract lines allowing you to also search by your part number for easy and accurate ordering of what you need. Your part number will appear in red in brackets directly after the Rowe code.

FF2900 (10617517) F
(MN1672) (472032)
Unit:100 Local Stock:64
Est. Lead days if no local or

See your agreed pricing

You will see your agreed contract and special pricing. This is indicated by the words **Your contract/special price** above the price.

Your Special price
\$10.00
VOLUME PRICING
QTY

Availability

You can check the availability and estimated lead times on products. Faster lead times may be available on request. The national stock is the total of our stock across the country.

TRACEABLE, IN CONTAMINATION PRO (R608+)
Unit:500mL Local Stock:463 National Stock:76
Est. Lead days if no local or national stock:5
NIS C.O.F.A. ADD TO QUICKLIST

Account manager use.

USERNAME: _____ PASSWORD: _____

If you don't have a Rowe Scientific login, go to www.rowe.com.au and click the register tab to request one.



www.rowe.com.au

Online 24 hours 7 days a week, by phone or face to face,
we give you the choice.



South Australia & NT
Ph: (08) 8186 0523
rowesa@rowe.com.au

Queensland
Ph: (07) 3376 9411
roweqld@rowe.com.au

Victoria & Tasmania
Ph: (03) 9701 7077
rowevic@rowe.com.au

Western Australia
Ph: (08) 9302 1911
rowewa@rowe.com.au

New South Wales & ACT
Ph: (02) 9603 1205
rowensw@rowe.com.au



*A 100% Australian
owned company
since 1987.*



REF-547 Prices do not include GST. We reserve the right to change specifications, details and descriptions without notice. Pictures for illustrative purposes only. Discounts do not apply to service, freight and or repairs. X\MARKETING\Labikes\JUL-SEPI\Labikes_Issues\OCT\DEC-19